

*Stand safely,
turn safely*



Etac Turner PRO

Safe transfers that are kind to
a carers back

etac[®]

Etac Turner PRO

Patient turning with standing support

The importance of safe transfers

People in need of assistance often feel unsafe in transfer situations, and caregivers are subject to strain when assisting with the transfer. The risk of accidents and insecurity from both caregiver and the assisted person may result in fewer transfers.

The Turner PRO allows both the assisted person and the caregiver to feel safe and the transfer performed with ease.

Easy to use

Turner PRO is the perfect tool for hospitals, care homes and in the home. The intuitive design makes it easy to use, the only adjustments needed are handle height and height of leg support.

The Turner PRO works well on different flooring materials including carpets.



A safe transfer which allows eye contact and good communication.

Turner PRO

Code 16090106

Dimensions

Handle height: 78-120 cm / 30.7"-49.2"

Footplate diameter: 40 cm / 15.7"

Footplate height: 20 mm / 0.8"

Lower leg support width: 38 cm / 15"

Lower leg support height: 15 cm / 5.9"

Material

Cast aluminum, powder coated steel, chrome 3

Compliant with

MDD 93/42/EEC

EN 12182

Total weight

8 kg / 17.6 lbs

Weight

Max user weight 200 kg / 441 lbs



Clean the product with a solvent-free cleaning agent (pH level 5–9), or with a 70% disinfectant solution. Can be decontaminated at max 85°C / 185°F if so required.




Features



Compact

Easy to transport. Easy to grip. Light, evenly balanced weight.

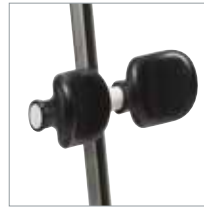

Angled leg support

Soft and comfortable, provides lateral support. Easy to adjust to whatever height is needed.



Contrasting coloured handle

The contrasting colour give the handle extra visibility. Offers many grip possibilities for both carer and user.



Functional base plate


Cast aluminum base plate with integrated anti slip function. Very easy to turn. The wheels make Turner PRO easy to transport indoors.



Design, function and safety are extremely important to us, and this is presented here in the form of a sit to stand aid. The scandinavian design gives you maximum comfort and a more stylish environment.



Turner PRO has a small footprint which makes it the perfect choice when transferring in a tight space, for example between a wheelchair and toilet.

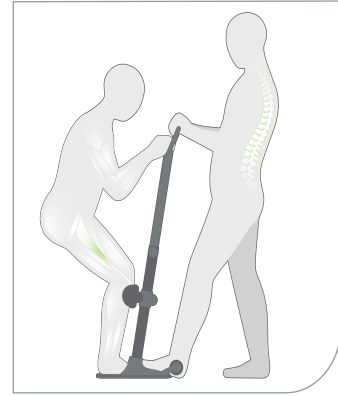

Encourages natural standing and sitting movement for the user. Enables a good posture and minimises the risk of injury for the carer when transferring.



Intended use

Turntable with support handle

- Transfer from bed to wheelchair / shower chair
- Transfer from wheelchair to toilet
- Transfer from chair to wheelchair / shower chair



*"Everyone should be able to live a free and independent life
and pursue their dreams regardless of any physical circumstances."*

*For the latest news and continuously
updated product information – please
visit **www.etac.com***

Distributed by **novis** MAY 2017 ©

Pressure care and patient handling specialists
novis.com.au / 1300 738 885